

# Miscellaneous

# Laser Glasses

TEKNOMEK®

www.teknomek.net




- Laser Glasses is a disposable product.
- It has %99 UV protection
- It is designed to protect eye from light, wind or any external impacts
- Com24 hours continuous use due to its flexibility.
- Logos and/or names of the hospitals can be printed if desired on to the film and/or elastic rubber



**AS101** Universal Eye Shield

Clear lightweight plastic. Features with smooth edges and ventilation holes for patient comfort. For use on left or right eye. Packaged sterile.

10 Per box

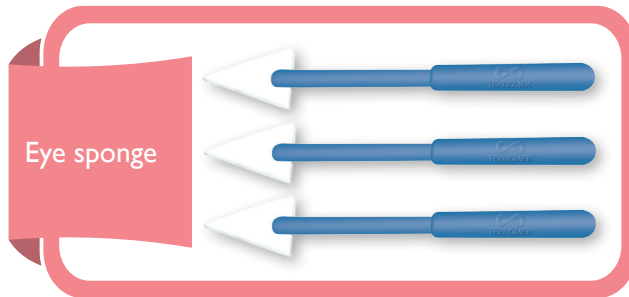


**AS200** Eye sponge (5 per pouch)

**AS201** Eye sponge (10 per pouch)

Made from highly absorbent, natural, cellulose material. Fast wicking action delivers quick absorption of fluid from the surgical area. The flat tip provides precise application at the surgical incision. 7cm overall length.


10 Per box

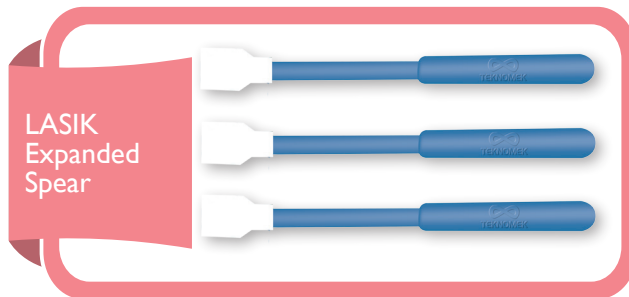

**AS210** LASIK Expanded Spear

Made from highly absorbent PVA material. Pre-expanded for convenience. Designed to provide instantaneous wicking. No need to pre-soak and fan.


10 Per box

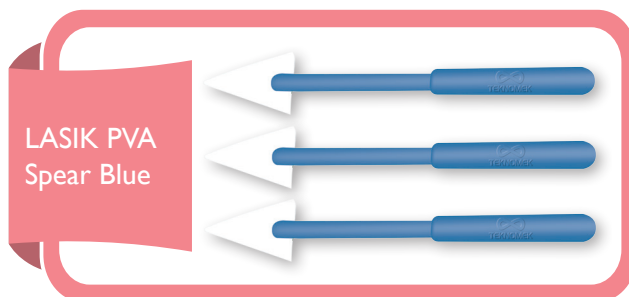

**AS211** LASIK PVA Spear Blue

Made from highly absorbent PVA material. Designed for optimal tissue manipulation and controlled fluid absorption to facilitate flap alignment and adhesion. Does not leave any particles.

10 Per box

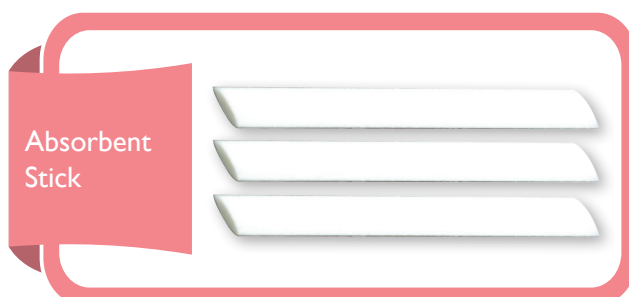

**AS220** Absorbent Stick 5 per pouch

**AS221** Absorbent Stick 10 per pouch

6.6cm overall length x 5mm diameter.

10 Per box

Cotton Swab



**AS230 Cotton Swab**

(Rounded/Pointed Tips), 7.6cm overall length. Rounded end is made for scleral depression and pointed end may be used as a second instrument to move sutures, manipulate tissue or absorb blood and fluid. Tightly wound ends create less linting than other cotton swabs.

10 Per box

High Flow Filter



**AS300 High Flow Filter**

Male/female luer lock fittings. Place between syringe and needle or cannula. Filters out particulates and microbiological contamination of .22 micron and larger.

10 Per box

High Flow Filter



**AS301 High Flow Filter with Syringe**

10 Per box

Instrument Wipe



**AS350 Instrument Wipe**

**AS351 Instrument Wipe**

**AS352 Instrument Wipe**

Wipes are designed to provide the surgeon with superior fiber free products for unparalleled fluid control and safe cleaning of fine microsurgical instruments. The exclusive material is highly absorbent, fast wicking, lint free, polyvinyl acetal (PVA). Available desired size.

10 Per box


Superflex Extension Set



**AS360 Superflex Extension Set**


25cm soft flexible tubing with male/female luer connectors.

10 Per box

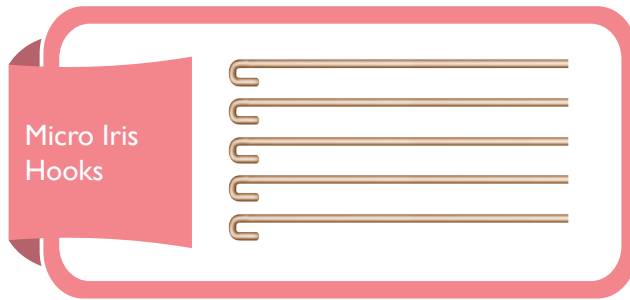


**AS370** Micro Iris Hooks

8mm long x .08mm diameter. Flexible, gold wire hooks designed to retract the iris during surgical procedures.  
5 Per box

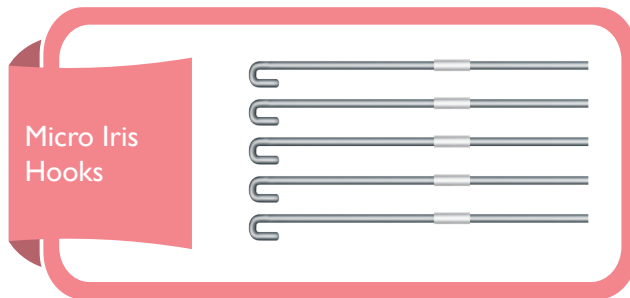


**AS371** Micro Iris Hooks

8mm long x .08mm diameter. Flexible wire hooks with silicone stoppers.  
5 Per box

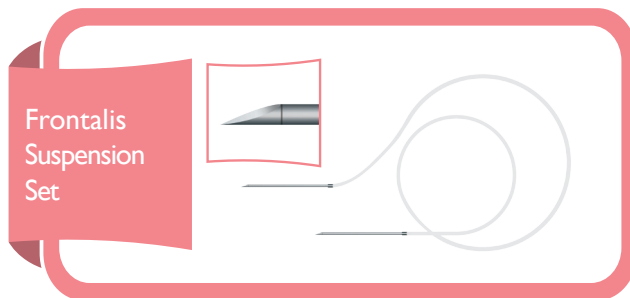

**AS380** Frontalis Suspension Set

Stainless steel needles .90mm x 6.3cm. Solid silicone tubing length 40cm. Silicone tubing sleeve length 10mm  
5 Per box

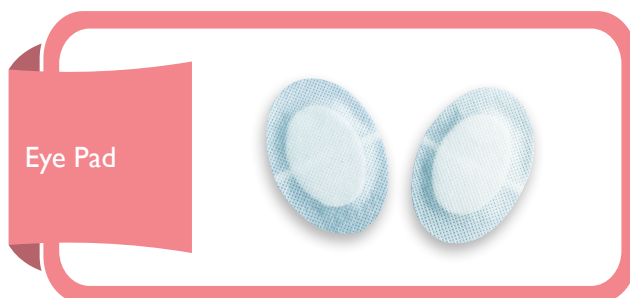


**AS390** Eye Pad

10 Per box

## Speculums

Eye Speculum (Adult Size)

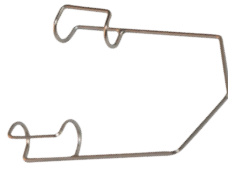


**AS400** Eye Speculum (Adult Size)

Eye Speculum [Barraquer] Adult, closed loop wire 14mm blades. Width 28mm. 0.8mm wire size.

10 Per box

Eye Speculum (Adult Size)

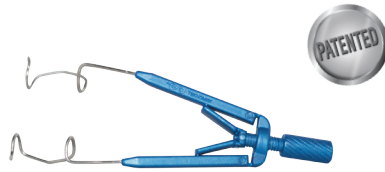


**AS400AI** Eye Speculum® (Adult Size)

Overall size 44 x 24mm, retraction size 14 x 5mm.

10 Per box

## Pediatric Speculums

Pediatric Eye Speculum

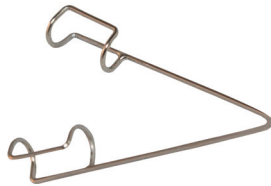

**AS410** Pediatric Eye Speculum (Small Size)

Overall size 30 x 21mm, retraction size 4 x 4mm.

**AS411** Pediatric Eye Speculum (Medium Size)

Overall size 34 x 22mm, retraction size 4.5 x 4.5mm.

**AS412** Pediatric Eye Speculum (Large Size)

Overall size 35 x 23mm, retraction size 10 x 5mm.

Pediatric Eye Speculum



**AS410AI** Pediatric Eye Speculum® (Small Size)

Overall size 30 x 21mm, retraction size 4 x 4mm.

**AS411AI** Pediatric Eye Speculum® (Medium Size)


Overall size 34 x 22mm, retraction size 4.5 x 4.5mm.

**AS412AI** Pediatric Eye Speculum® (Large Size)

Overall size 35 x 23mm, retraction size 10 x 5mm.

## ROP Sets

ROP Set

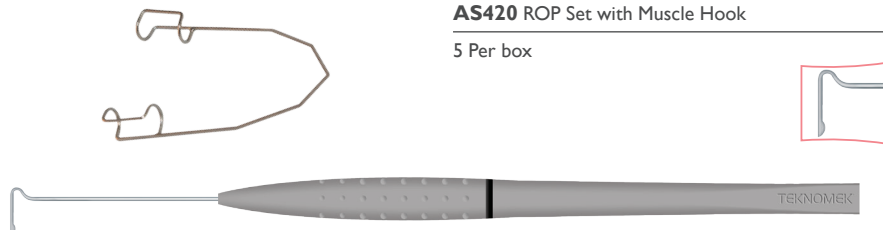


**AS420** ROP Set with Muscle Hook

5 Per box



**ROP Set**


**AS420P ROP Set with Muscle Hook**  
5 Per box

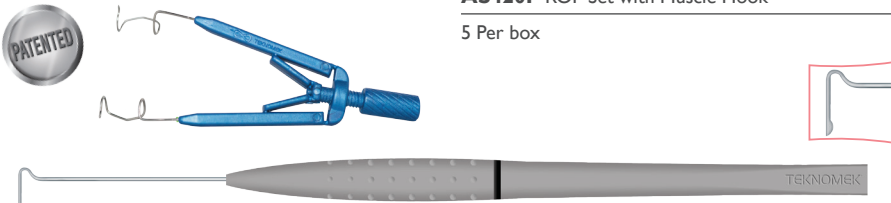

The image shows a blue muscle hook and a grey handle with a silver tip. A circular 'PATENTED' logo is visible. A small inset image shows the hook's profile.

**ROP Set with T Hooks**


**AS430 ROP Set with T Hooks**  
5 Per box

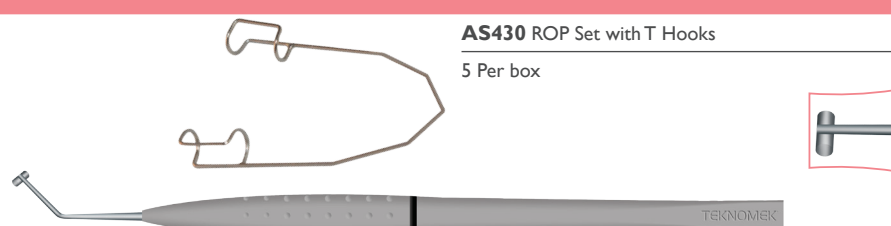

The image shows a grey handle with a silver tip and a wire T-hook. A small inset image shows the T-hook's profile.

**ROP Set with T Hooks**


**AS430A ROP Set with T Hooks**  
5 Per box

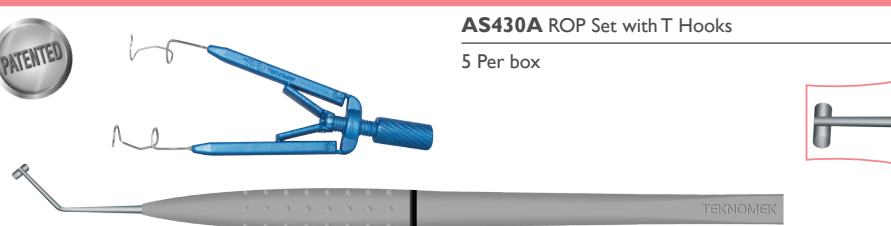

The image shows a blue muscle hook and a grey handle with a silver tip. A circular 'PATENTED' logo is visible. A small inset image shows the T-hook's profile.

**ROP Set with Ans Hooks**


**AS440 ROP Set with Ans Hooks**  
5 Per box

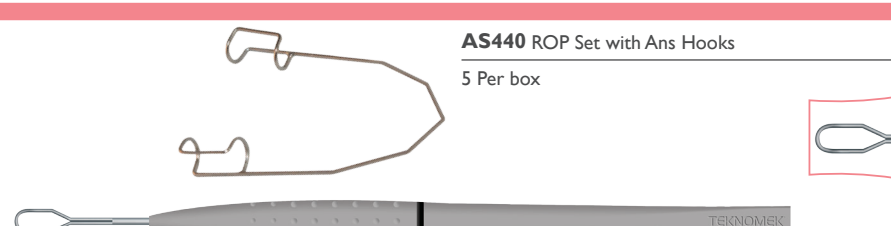

The image shows a grey handle with a silver tip and a wire Ans hook. A small inset image shows the Ans hook's profile.

**ROP Set with Ans Hooks**

**AS440A ROP Set with Ans Hooks**  
5 Per box

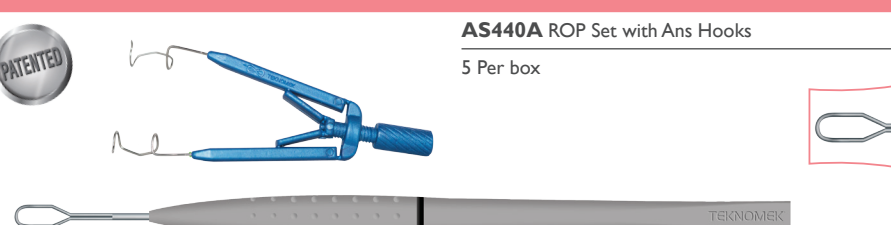

The image shows a blue muscle hook and a grey handle with a silver tip. A circular 'PATENTED' logo is visible. A small inset image shows the Ans hook's profile.

## Riboflavin



### Riboflavin


- Isotonic ophthalmic solution which contains %0.1 Riboflavin
- Used in combination with UV-A light after removal of corneal epithelium which enables collagen cross-linking
- Protects deeper tissue penetration (corneal endothelium, lens and retina) of radiation (up to 95%) acting as a barrier
- Used especially in documented progressive keratoconus and also for other corneal ectatic disorders

### Riboflavin TE

- Hypotonic ophthalmic solution which contains %0.1 Riboflavin
- Used in combination with UV-A light over the intact epithelium which enables collagen cross-linking
- Transepithelial penetration aided by the addition of enhancers (eliminating the need for corneal epithelial removal) reduces procedure time as well as postoperative discomfort
- Protects deeper tissue penetration (corneal endothelium, lens and retina) of radiation (up to 95%) acting as a barrier
- Used especially in documented progressive keratoconus and also for other corneal ectatic disorders

## Product Properties

	Riboflavin Hypotonic	Riboflavin TE	Riboflavin
Basic Riboflavin	0,100 g	0,100 g	0,100 g
Dextran T500	-	15 g	20 g
Purified Water	q.b. a 100 g	q.b. a 100 g	q.b. a 100 g

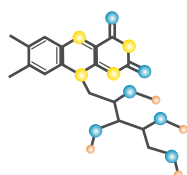

ml	Riboflavin Hypotonic	Riboflavin TE	Riboflavin
0.5ml	RBH05	RBSTE05	RBS05
1ml	RBH1	RBSTE1	RBS1
2ml	RBH2	RBSTE2	RBS2
3ml	RBH3	RBSTE3	RBS3

Product Description : Presented in Syringe

Avoid from light exposure

Store at room temperature